

# MAPKAPK-2 (35-I): sc-100393

## BACKGROUND

The MAPKAP kinases (for MAP kinase activated protein kinases) are a group of MAP kinase substrates which are themselves kinases. In response to activation, the MAP kinases phosphorylate downstream components on a consensus Pro-X-Ser/Thr-Pro motif. Several kinases that contain this motif have been identified and serve as substrates for the ERK and p38 MAP kinases. These include the serine/threonine kinases Rsk-1 (also designated MAPKAP kinase-1), Rsk-2 and Rsk-3, which are phosphorylated by ERK 1 and ERK 2. Similarly, p38 phosphorylates and activates the serine/threonine kinases MAPKAP kinase-2 and MAPKAP kinase-3 (also designated 3pK). The serine/threonine kinases Mnk1 and Mnk2 are substrates for both ERK and p38 MAP kinases.

## REFERENCES

1. Sturgill, T.W., et al. 1988. Insulin-stimulated MAP2 kinase phosphorylates and activates Ribosomal Protein S6 kinase II. *Nature* 334: 715-718.
2. Stokoe, D., et al. 1992. MAPKAP kinase-2: a novel protein kinase activated by mitogen-activated protein kinase. *EMBO J.* 11: 3985-3994.
3. Davis, R.J. 1993. The mitogen-activated protein kinase signal transduction pathway. *J. Biol. Chem.* 268: 14553-14556.

## CHROMOSOMAL LOCATION

Genetic locus: MAPKAPK2 (human) mapping to 1q32.1; Mapkapk2 (mouse) mapping to 1 E4.

## SOURCE

MAPKAPK-2 (35-I) is a mouse monoclonal antibody raised against recombinant MAPKAPK-2 of human origin.

## PRODUCT

Each vial contains 100 µg IgG<sub>2b</sub> kappa light chain in 1.0 ml of PBS with < 0.1% sodium azide and 0.1% gelatin.

## APPLICATIONS

MAPKAPK-2 (35-I) is recommended for detection of MAPKAPK-2 of mouse, rat and human origin by Western Blotting (starting dilution 1:200, dilution range 1:100-1:1000), immunoprecipitation [1-2 µg per 100-500 µg of total protein (1 ml of cell lysate)], immunofluorescence (starting dilution 1:50, dilution range 1:50-1:500), immunohistochemistry (including paraffin-embedded sections) (starting dilution 1:50, dilution range 1:50-1:500) and solid phase ELISA (starting dilution 1:30, dilution range 1:30-1:3000).

Suitable for use as control antibody for MAPKAPK-2 siRNA (h): sc-35855, MAPKAPK-2 siRNA (m): sc-35856, MAPKAPK-2 shRNA Plasmid (h): sc-35855-SH, MAPKAPK-2 shRNA Plasmid (m): sc-35856-SH, MAPKAPK-2 shRNA (h) Lentiviral Particles: sc-35855-V and MAPKAPK-2 shRNA (m) Lentiviral Particles: sc-35856-V.


Molecular Weight of MAPKAPK-2: 45 kDa.

Positive Controls: A-431 whole cell lysate: sc-2201, K-562 whole cell lysate: sc-2203 or NIH/3T3 whole cell lysate: sc-2210.


## RECOMMENDED SUPPORT REAGENTS

To ensure optimal results, the following support reagents are recommended: 1) Western Blotting: use m-IgGκ BP-HRP: sc-516102 or m-IgGκ BP-HRP (Cruz Marker): sc-516102-CM (dilution range: 1:1000-1:10000), Cruz Marker™ Molecular Weight Standards: sc-2035, UltraCruz® Blocking Reagent: sc-516214 and Western Blotting Luminol Reagent: sc-2048. 2) Immunoprecipitation: use Protein A/G PLUS-Agarose: sc-2003 (0.5 ml agarose/2.0 ml). 3) Immunofluorescence: use m-IgGκ BP-FITC: sc-516140 or m-IgGκ BP-PE: sc-516141 (dilution range: 1:50-1:200) with UltraCruz® Mounting Medium: sc-24941 or UltraCruz® Hard-set Mounting Medium: sc-359850. 4) Immunohistochemistry: use m-IgGκ BP-HRP: sc-516102 with DAB, 50X: sc-24982 and Immunohistomount: sc-45086, or Organo/Limonene Mount: sc-45087.

## DATA

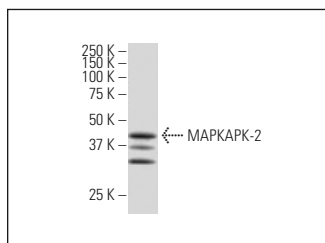

MAPKAPK-2 (35-I): sc-100393. Western blot analysis of MAPKAPK-2 expression in K-562 whole cell lysate.

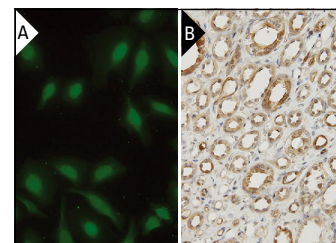

MAPKAPK-2 (35-I): sc-100393. Immunofluorescence staining of paraformaldehyde-fixed HeLa cells showing nuclear localization (A). Immunoperoxidase staining of formalin-fixed, paraffin-embedded human kidney tissue showing cytoplasmic localization (B).

## SELECT PRODUCT CITATIONS

1. Garbati, M.R., et al. 2013. FANCA and FANCC modulate TLR and p38 MAPK-dependent expression of IL-1β in macrophages. *Blood* 122: 3197-3205.

## STORAGE

Store at 4° C, \*\*DO NOT FREEZE\*\*. Stable for one year from the date of shipment. Non-hazardous. No MSDS required.

## RESEARCH USE

For research use only, not for use in diagnostic procedures.

## PROTOCOLS

See our web site at [www.scbt.com](http://www.scbt.com) for detailed protocols and support products.